


QIMAGING QUICK INSTALLATION GUIDE


Install the Camera-Interface Software

1. Install the high speed PCIe card in your computer.
2. Power on the computer and insert the QImaging USB stick into your USB drive.
3. Install camera drivers by choosing "QI_setup.exe".
4. Follow the installation instructions on screen.
5. Restart your computer and connect the camera to the PC.

optiMOS™

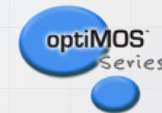

1. Plug one end of the data cable into a high speed PCIe port of your computer. Plug the other end of the cable into the camera's data socket.
2. Connect the power cord to the power supply, plug the power cord into a wall outlet, and connect the power supply to the camera.
3. Power on your computer. Then power on the camera by pushing the power button.
NOTE: "WAIT" LED should be off before running the camera. It takes approximately 30 seconds for the camera to initialize.

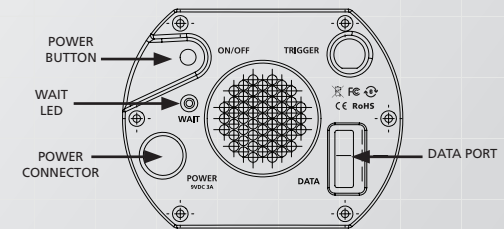

More detailed installation instruction can be found in the camera instruction manual on the USB stick.

Contact QImaging Customer Service for additional information: 1-800-874-9789 or www.qimaging.com



www.qimaging.com

©2015 QImaging. All rights reserved.

58-607-001 Rev A0